

# Crystal Digital PCR® Assay

## Information Sheet

For Research Use Only. Not for use in diagnostic procedures.

### Product Name

MLH1 Methylation Crystal Digital PCR® Assay (R51036)

### Description

#### Detected Targets

Targets	Sample Type	Detection Channels	Multiplex
<b>MLH1 Methylated/Unmethylated Promoter</b>	DNA	Blue/Green/Red	3

The MLH1 Methylation Crystal Digital PCR® Assay is a 10X assay designed to detect and quantify methylation in the promoter of MLH1 gene using the Ruby Chip. MLH1 is involved in DNA mismatch repair and promoter hypermethylation is the most common reason of gene inactivation. This is an important cause for genomic instability and cell proliferation in several cancers, notably in colorectal and endometrial cancers.

#### Multiplexing Strategy: Color-Combination

This assay relies on the Color-Combination multiplexing strategy proprietary to Stilla Technologies, in which targets are detected, differentiated, and quantified by Crystal Digital PCR® using 2 fluorophores.

The table below indicates with a “X” which channel(s) are used for each target in the assay:

Targets	Blue	Teal	Green	Yellow	Red	Infra-Red	Long-Shift
<b>MLH1 Unmethylated-DNA</b>	X				X		
<b>MLH1 Methylated-DNA</b>	X		X				

### Components

MLH1 Methylation Crystal Digital PCR® Assay comprises two reagents: a pool of the assay specific primers and Crystal Flex Probes and a Positive Control. Please refer to the lot specific Certificate of Conformity for characterized concentration, available for download at the Technical Resources section of the Stilla Technologies website.

Component Name	Reference	Concentration	Description
<b>MLH1 Methylation Crystal Digital PCR® Assay</b>	R51036	10X	Detects methylation of the MLH1 promoter
<b>MLH1 Positive Control</b>	R51036.PC0	10X	Contains: synthetic sequences corresponding to methylated and unmethylated DNAs after bisulfite treatment

## Specific Recommendation Regarding Sample Treatment and DNA Input

The assay is designed to detect methylated and non-methylated sequences after bisulfite treatment. Samples must therefore first be subjected to bisulfite treatment. The kit used during assay validations is indicated in section “Consumables Required but Not Provided”.

## Thermocycling Programs

### On the naica® system:

	Step	Ramp rate
<b>Step 1</b>	Partition for Ruby Chip	-
<b>Step 2</b>	Temperature 95°C for 180 seconds	1°C/sec
<b>Step 3</b>	Begin Loop for 60 Iterations	-
<b>Step 3.1</b>	Temperature 95°C for 15 seconds	1°C/sec
<b>Step 3.2</b>	Temperature 60°C for 60 seconds	1°C/sec
<b>Step 4</b>	Temperature 58°C for 300 seconds	1°C/sec
<b>Step 5</b>	Release for Ruby Chip	-

### On the Nio™ Digital PCR:

	Step	Ramp rate
<b>Step 1</b>	Partition for Ruby Chip	-
<b>Step 2</b>	Temperature 95°C for 180 seconds	1°C/sec
<b>Step 3</b>	Begin Loop for 60 Iterations	-
<b>Step 3.1</b>	Temperature 95°C for 15 seconds	2°C/sec
<b>Step 3.2</b>	Temperature 60°C for 60 seconds	2°C/sec
<b>Step 4</b>	Temperature 58°C for 300 seconds	1°C/sec
<b>Step 5</b>	Release for Ruby Chip	-

## Image Acquisition

Download the dedicated scanning file from the Technical Resources section of the Stilla Technologies website:

- ScanningTemplate\_Prism3\_MLH1\_R51036.ncx (3-color naica® system)
- ScanningTemplate\_Prism6\_MLH1\_R51036.ncx (6-color naica® system)
- NioProtocol\_3C-60X-60°C-60s.nioprotocol (Nio™ Digital PCR)
- NioAssay\_3C\_MLH1\_R51036.nioassay (Nio™ Digital PCR)

## Image Analysis

The following files are embedded in the dedicated scanning files listed above:

- CompensationMatrix\_Prism3\_MLH1\_R51036.ncm (3-color naica® system)
- UniversalCompMatrix\_3C\_Prism6-Nio.ncm (6-color naica® system, Nio™ Digital PCR)
- AnalysisConfiguration\_MLH1\_R51036.nca (all systems)

## Consumables Required but Not Provided

- Ruby Chip (C16011)
- naica® PCR MIX 10X (R10106)
- Crystal Universal Reporters 3 (R41401 200 reactions, R41402 1000 reactions)
- Nuclease-free water
- Bisulfite conversion kit (Example: EZ DNA Methylation-Gold Kit, ref: ZD5005 or ZD5006 from Ozyme)

## Instruction for PCR Mix Preparation

Specific instructions for preparing the PCR mix are given below.

Reagent Name		Initial Concentration	Final Concentration	Volume per reaction (µL)
naica® PCR MIX Buffer A	●	10x	1x	0.60
naica® PCR MIX Buffer B	●	100%	4%	0.24
Crystal Digital PCR® Assay	●	10x	1x	0.60
Crystal Universal Reporter Tube A	●	40x	1x	0.15
Nuclease-free water		NA	NA	Variable
<b>Template DNA</b>		<b>NA</b>	<b>NA</b>	<b>Variable</b>
<i>or Positive Control</i>	○	10x	1x	0.60
<b>Total reaction volume (µL)</b>				<b>6.0</b>

## Representative Data and Instructions for Analysis

Set thresholds for separating positive and negative populations on the 1D plots. To optimize the analysis, the Blue/Green/Red thresholds should be set above the negative cluster. Examples of results obtained on the Nio™+ system are given below.

Remark: The threshold can be adjusted on each individual chamber to optimize its placement.

Wet lab testing was carried out using H<sub>2</sub>O as a negative control and a positive control consisting of synthetic DNAs corresponding to methylated and unmethylated DNA sequences after bisulfite treatment. Synthetic DNAs, hgDNA and methylated DNA standard (CpG Methylated Human Genomic DNA, ref: SD1131 from ThermoFisher) were also tested individually after bisulfite treatment.

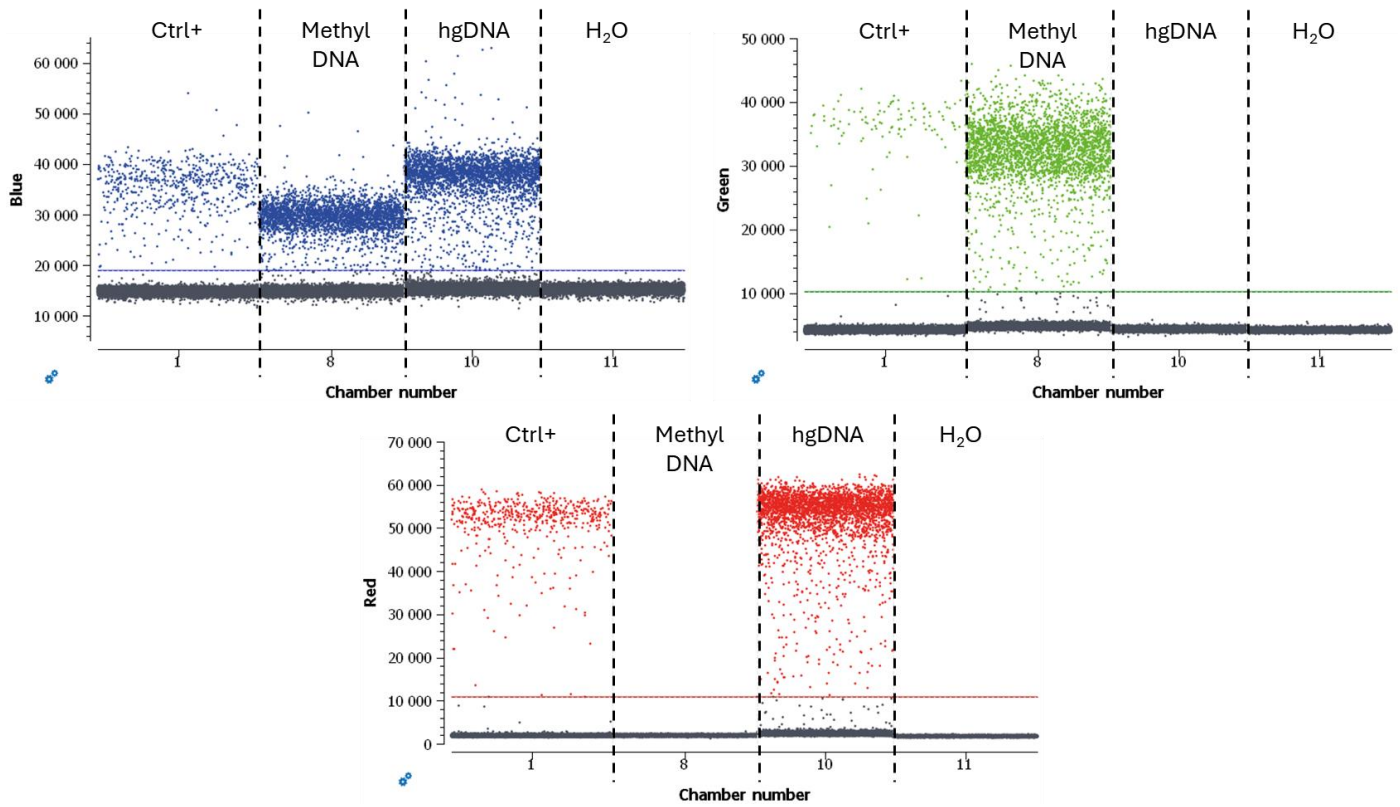


Figure 1: 1D plots obtained during wet lab testing on the Nio™+. The thresholds should be set above negative cluster.



Stilla Technologies  
F-94800 Villejuif, FRANCE

Registered names and trademarks used in this document, even when not specifically marked, are not to be considered unprotected by law.